

AERObag®

Silicone and PVC Manual Resuscitators



Manual Resuscitators

Resuscitation Masks

PEEP valves

Accessories

Reservoirvalve integrated - "All-In-One"

The next generation



Powerful helpers for patients and user friendly resuscitation

Flexibility

3 models (bag volumes): adult (1500 ml) | children (550 ml) | infant (280 ml).

Different masks and sizes available to accommodate a wide range of patients needs

Silicone • PVC

3 Models each

Models

Silicone-Manual Resuscitators/

Adults | Children | Infants



Models

PVC-Manual Resuscitators

Adults | Children | Infants



High-quality materials

Long lifespan! Operational even in toughest emergency conditions.

- Autoclavable up to 134°C (Manual Resuscitators of silicone and polysulfone)
- Easy to disinfect and sterilize by gas (Manual Resuscitators of PVC and polycarbonate)

High ease of use

- Comfortable “grip”
- Unique right-angle connector swivels in any direction
- Pressure relief valve 60 mbar - adults (fig. 1) or 40 mbar - children and infants
- Can be folded to save space

Hygienic and transparent

- Contamination easy to detect
- Easy to clean

Exact fit

- Connections according to ISO 5356
- Masks, tubings and filters can be connected
- Innovative All-In-One inlet-valve (fig. 2)
 - ➔ integrated Reservoir valve ①
- Oxygen gas inlet for 6 mm tubing ②
- Direct connection for Reservoir bag, demand valve respectively bag-refill systems (conic plug - 26mm OD / 22 mm ID) ③



Abb. 1

Inlet valve (integrated reservoir valve)

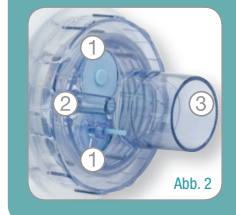


Abb. 2

AERObag® - Resuscitation Masks

Anatomically shaped resuscitation masks seal well at low pressure.

Different mask sizes for adults, children and infants.

Silicone Resuscitation Masks (reusable / autoclavable)

- Masks with silicone cuff and a clear rigid shell of polysulfone ④
- Full-silicone masks with integrated cuff ⑤

PVC Resuscitation Masks (single use/ disposable)

- Masks with inflatable cuff (with valve) and clear rigid shell of PVC
→ colour coded fixing hook of PE (for fixing strip) ⑥
- Full-PVC masks with integrated cuff ⑦

4 different types of masks




AERObag® - PEEP-Valves (PEEP = Positive Endexpiratory Pressure)

Improvement in oxygenation at varied disturbances of pulmonary function.

All PEEP valves are equipped with 2 adapters with the dimensions 22 mm OD and 30 mm ID ⑧.

They will be attached to the patient valve via diverter/adapter (Silicone ⑨ a | PVC ⑨ b).

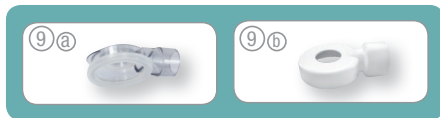

PEEP-Valves

Reusable / autoclavable ⑧ a

- PEEP valve 10 mbar - continuously adjustable 2 cm H₂O - 10 cm H₂O
- PEEP valve 10 mbar - continuously adjustable 5 cm H₂O - 20 cm H₂O

Single use / disposable ⑧ b

- PEEP valve 10 mbar - continuously adjustable 2 cm H₂O - 10 cm H₂O
- PEEP valve 10 mbar - continuously adjustable 5 cm H₂O - 20 cm H₂O

2 Models each



AERObag® - Storage Box

Everything stored and quick to release.

Storage Box suitable for all AERObag® models incl. accessories.

- Clear plastic
- Hinged handle
- Large storage space



AERObag® - Silicone - Manual Resuscitators

Item-No.	Description
HBB 05-E	AERObag® - BB05-E, Manual Resuscitator, silicone, adults
HBB 05-K	AERObag® - BB05-K, Manual Resuscitator, silicone, children
HBB 05-B	AERObag® - BB05-B, Manual Resuscitator, silicone, infants



Spare parts

Item-No.	Description
HBB 05-PVD40	AERObag® - PVD05-40, patient valve with pressure relief valve 40 cmH ₂ O, complete
HBB 05-PVD60	AERObag® - PVD05-60, patient valve with pressure relief valve 60 cmH ₂ O, complete
HBB 05-PVD40-DB	pressure relief valve 40 cmH ₂ O, single
HBB 05-PVD60-DB	pressure relief valve 60 cmH ₂ O, single
HBB 05-PVD-OT	upper part without pressure relieve valve
HBB 05-PVD-LM	membrane, blue
HBB 05-PVD-UT	bottom part
HBB 05-PVD-AAM	exhalation membrane, blue
HBB 05-HR	retainer sealing ring

Balg AERObag® - BE05, BK05 und BB05

Item-No.	Description
HBB 05-BE	AERObag® - BE05, bag adults, without valves, autoclavable
HBB 05-BK	AERObag® - BK05, bag children, without valves, autoclavable
HBB 05-BB	AERObag® - BB05, bag infants, without valves, autoclavable

Ansaugventil AERObag® - AVRV05

Best.-Nr.	Description
HBB 05-AVRV	AERObag® - AVRV05, inlet valve with integrated reservoir valve threepart, complete, autoclavable

Zubehör

Item-No.	Description
HBB 05-RB600	AERObag® - RB06-600, reservoir bag 600 ml
HBB 05-RB2500	AERObag® - RB06-2500, reservoir bag 2500 ml
HBB 05-DV	AERObag® - DV05, Diverter / PEEP valve adapter
HSS 05-2	AEROpart® - O ₂ -tubing, 2 m
HBB 05-HS	AERObag® - HS05, hand strap, silicone



AERObag® - Silicone - Resuscitation Masks

Item-No.	Description
HBM 05-5	AERObag® - BM05-5, silicone resuscitation mask, sz. 5
HBM 05-4	AERObag® - BM05-4, silicone resuscitation mask, sz. 4
HBM 05-3	AERObag® - BM05-3, silicone resuscitation mask, sz. 3
HBM 05-5V	AERObag® - BM05-5V, full-silicone resuscitation mask, sz. 5
HBM 05-4V	AERObag® - BM05-4V, full-silicone resuscitation mask, sz. 4
HBM 05-3V	AERObag® - BM05-3V, full-silicone resuscitation mask, sz. 3
HBM 05-2V	AERObag® - BM05-2V, full-silicone resuscitation mask, sz. 2
HBM 05-1V	AERObag® - BM05-1V, full-silicone resuscitation mask, sz. 1
HBM 05-0V	AERObag® - BM05-0V, full-silicone resuscitation mask, sz. 0



AERObag® - PEEP Valves (reusable)

Item-No.	Description
HPV 05-10	AERObag® - PV05-10, PEEP valve, 10 mbar, autoclavable, with 2 adapter (22 mm AD; 30 mm ID), pressure continuously adjustable 2 cm - 10 cm H ₂ O
HPV 05-20	AERObag® - PV05-20, PEEP valve, 20 mbar, autoclavable, with 2 adapter (22 mm AD; 30 mm ID), pressure continuously adjustable 5 cm - 20 cm H ₂ O



AERObag® - PVC - Manual Resuscitators

Item-No.	Description
HBB 06-E	AERObag® - BB06-E, Manual Resuscitator, PVC, adults
HBB 06-K	AERObag® - BB06-K, Manual Resuscitator, PVC, children
HBB 06-B	AERObag® - BB06-B, Manual Resuscitator, PVC, infants



Spare parts

Patientenventil AERObag® - PVD06

Item-No.	Description
HBB 06-PVD40	AERObag® - PVD06-40, patient valve with pressure relief valve 40 cmH ₂ O, complete
HBB 06-PVD60	AERObag® - PVD06-60, patient valve with pressure relief valve 60 cmH ₂ O, complete
HBB 06-PVD40-DV	AERObag® - PVD06-40DV, patient valve with pressure relief valve 40 cmH ₂ O, complete, with diverter / port 30 mm AD) for PEEP valves
HBB 06-PVD60-DV	AERObag® - PVD06-60DV, patient valve with pressure relief valve 60 cmH ₂ O, complete, with diverter / port 30 mm AD) for PEEP valves

Balg AERObag® - BEK06, BKK06 und BBK06:

Item-No.	Description
HBB 06-BEK	AERObag® - BE06, bag adults, with inlet valve and retainer
HBB 06-BKK	AERObag® - BK06, bag children, with inlet valve and retainer
HBB 06-BBK	AERObag® - BB06, bag infants, with inlet valve and retainer

Zubehör

Item-No.	Description
HBB 06-RB600	AERObag® - RB06-600, reservoir bag 600 ml
HBB 06-RB2500	AERObag® - RB06-2500, reservoir bag 2500 ml
HBB 06-DV30	AERObag® - DV06-30, Diverter / PEEP valve adapter
HSS 05-2	AEROpart® - O ₂ -tubing, 2 m
HBB 06-HS	AERObag® - HS06, hand strap, PVC



AERObag® - PVC - Resuscitation Masks

Item-No.	Description
HBM 06-5V	AERObag® - BM06-5V, PVC resuscitation mask with valve, sz. 5
HBM 06-4V	AERObag® - BM06-4V, PVC resuscitation mask with valve, sz. 4
HBM 06-3V	AERObag® - BM06-3V, PVC resuscitation mask w. valve, sz. 3
HBM 06-2V	AERObag® - BM06-2V, PVC resuscitation mask w. valve, sz. 2
HBM 06-1V	AERObag® - BM06-1V, PVC resuscitation mask w. valve, sz. 1
HBM 06-0V	AERObag® - BM06-0V, PVC resuscitation mask w. valve, sz. 0
HBM 06-5E	AERObag® - BM06-5E, full-PVC resuscitation mask, sz. 5
HBM 06-4E	AERObag® - BM06-4E, full-PVC resuscitation mask, sz. 4
HBM 06-3E	AERObag® - BM06-3E, full-PVC resuscitation mask, sz. 3



AERObag® - PEEP Valves (single use/disposable)

Item-No.	Description
HPV 06-10	AERObag® - PV06-10, PEEP valve, 10 mbar, single use, with 2 adapter (22 mm AD; 30 mm ID), pressure continuously adjustable 2 cm - 10 cm H ₂ O
HPV 06-20	AERObag® - PV06-20, PEEP valve, 20 mbar, single use, with 2 adapter (22 mm AD; 30 mm ID), pressure continuously adjustable 5 cm - 20 cm H ₂ O



Aufbewahrungsbox AERObag® - BB05-AB

Item-No.	Description
HBB 05-AB	AERObag® - BB05-AB, storage box for manual resuscitators plastic, internal dimension (W,H,D) 310 x 110 x 130 mm

HUM

HUM Gesellschaft für Homecare
und Medizintechnik mbH

Zum Pier 79 | D-44536 Lünen

FON +49 (0)2 31/ 88 08 85-0 | FAX +49 (0)2 31/ 88 08 85-58

Mail: sales@hum-online.de | http://www.hum-online.de



AERObag® - Manual Resuscitators

PP-AEROBAG_GB | Stand 21.10.2013 | Revision 2

© HUM Gesellschaft für Homecare und Medizintechnik mbH