

AEROSuc®

Foot Suction Pump | Hand-held Suction Pump



First Responder

Emergency Medical Service

Emergency Medicine

Medical Practices

Manual Suction Pumps

Powerful and helpful in every emergency situation

Foot Suction Pump

Quality at
its best

AEROSUC[®] - VacUp



Retaining clip with integrated vacuometer and for pressure reduction | ventilation

Dimensions (LxWxH): 180 x 130 x 110 mm

Clear advantages

Developed for daily use | no power required and free of maintenance

With 1300mm suction hose (Ø 8 mm) for suction of big accumulations of blood, saliva etc., and semifluid and solid parts of food

650ml suction canister with overfilling valve

Easy and certain handling

Can also be used for evacuating vacuum mattresses

Robust despite of little weight (850 g)

Powerful at maximum Pressure of -615 mmHg | -0,82 bar, with suction volume of 72 lpm. (Air)

High quality material | autoclavable at 121 °C

1 Retaining clip | 3 Functions

3 functions

Adjustment of the pedal
For space-saving storage

Vacuometer
For reading out the pressure on the scale

Venting device
For pressure reduction and venting by inserting the clip completely



Retaining clip



Hand-held Suction Pump

AEROsuc® - Rescue

Easy to
success




- Catheter adapter for funnel shape suction catheters
- Safe handling with low weight (275 g)
- Dimensions (HxWxD): 180 x 180 x 70 mm
fits in all current emergency bags | backpacks

- Clear advantages**
- Developed for the fast use; no power required and free of maintenance
 - With one suction tube size L (adult) and one suction tube size S (children)
 - 300 ml suction jar with overfilling valve
 - Easy one-hand use with maximum performance
 - Powerful at maximum Pressure of -550 mmHg | -0,73 bar, with suction volume of > 20 lpm. (Air)

Protection system devices

Catheter adapter

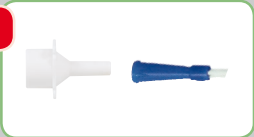

Overfilling valve

Contamination guard for the pump

Runback guard to the patient

Valve flap



Seal-system at pump




AEROSUC® - VacUp - Foot Suction Pump, complete

Art. Nr.: HAG 01-FP

Spare parts

Description	Material	Item-Number
① Vacuum pump	Polycarbonate (PC)	HAG01-FP-VP
② Filter	Tecnopor	HAG01-FP-F
③ Canister seal	Thermoplastic Elastomer (TPE)	HAG01-FP-D
④ Canister	Polycarbonate (PC)	HAG01-FP-AB
⑤ Ball	Polypropylene (PP)	HAG01-FP-K
⑥ Canister lid	Thermoplastic Elastomer (TPE)	HAG01-FP-AB-D
⑦ Tube	Polyvinylchloride (PVC)	HAG01-FP-S
Tubing connector	Polypropylene (PP)	HAG01-FP-SV

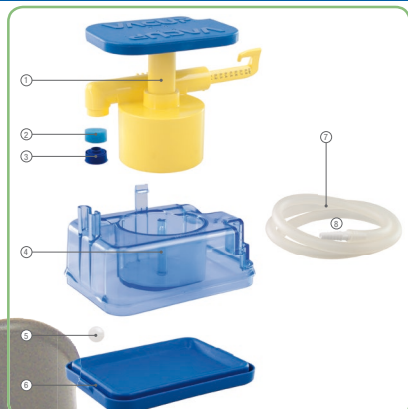

AEROSUC® - Rescue - Hand-held Suction Pump, complete

Art. Nr.: HAG 01-HP

Spare parts

Description	Item-Number
① Disposable suction jar, 300 ml	HAG01-HP-AG
② Disposable-PVC suction tubing, size L (adult)	HAG01-HP-ASL
③ Disposable-PVC suction tubing, size S (children)	HAG01-HP-ASS
④ Adapter for funnel shape suction catheter	HAG01-HP-AAK



HUM

HUM Gesellschaft für Homecare
und Medizintechnik mbH
Zum Pier 79 | D-44536 Lünen

FON +49 (0)2 31/ 88 08 85-0 | FAX +49 (0)2 31/ 88 08 85-58
Mail: sales@hum-online.de | <http://www.hum-online.de>

CE

AEROSUC® - Manual Suction Pumps
PP-AEROSUC-MAP_GB | Stand 05.11.2014 | Revision 1
© HUM Gesellschaft für Homecare und Medizintechnik mbH